

Noviomagus Mini Mill



Today's technology guarantees high quality chips

High quality construction

The Noviomagus Mini Mill is unique because of its superior design and manufacturing.


It only consists of few parts, a.o.: a housing, milling drums, a push block and a handle. All parts of the Noviomagus Mini Mill are made of machined stainless steel, which makes it an extremely accurate and durable device. There is no risk for abrasion. Additionally, the milling drums are made of a special hardened stainless steel alloy. Even in case of frequent use they do not show signs of wear. This combination guarantees a high quality graft without metal contamination.



Easy and optimal use

The teeth of the milling drums feature a unique patented technology. They are spirally positioned and razor sharp. Because of the stable base plate and optional clamp, the Noviomagus Mini Mill is easy to handle and it takes little strength to mill the bone parts. Because of the optimal design of the housing, all bone is milled without bone pieces remaining in the funnel. It can mill boneparts in less than 10 seconds.

spirally positioned

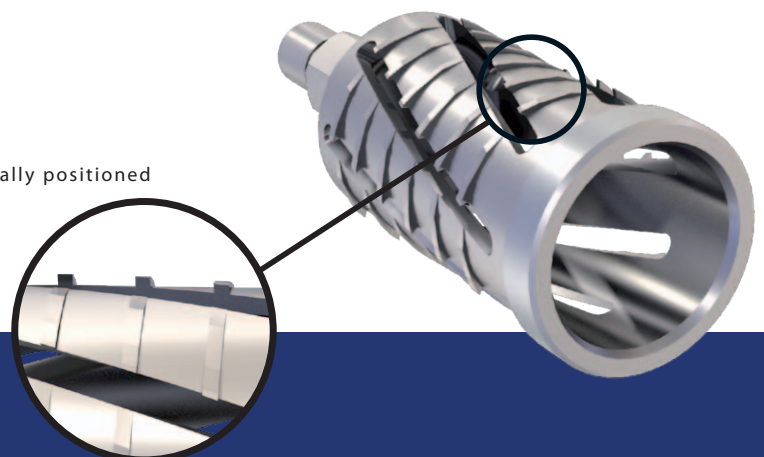

Can be used instantaneously

The Noviomagus Mini Mill guarantees bone chips with an intact trabecular bone structure, because it breaks out chips. The Noviomagus Mini Mill can be used for milling bone from the iliac crest, jaw and spine. The chips can be used instantaneously during surgery.

The Extra Fine milling drum can be used for spinal surgery. The Extra Extra Fine milling drum is designed to produce extra small bone chips which can be used for maxillo-facial and hand-foot surgery. The Extra Extra Extra Fine milling drum is meant for dental surgery, e.g. implantology.



Less than 10 seconds!



Advanced benefits

- Milling of bone particles takes less than 10 seconds
- 3 years guarantee for hospital use
- Three types of milling drums for small chip sizes
 - Extra Fine (3 - 5 mm)
 - Extra Extra Fine (1 - 3 mm)
 - Extra Extra Extra Fine (0.5 - 1 mm)
- No risk of metal abrasion
- Long life-span without wear
- Controlled bone particle size
- Easy to clean
- Optional clamp for additional stability



Today's technology guarantees high quality chips

Ordering information

Noviomagus Mini Mill


Noviomagus Mini Mill :

Mini Housing	SBM-1002-29
Mini Push Block	SBM-1002-30
Mini Handle	SBM-1002-31
Mini Screw Nut	SBM-1002-32
Noviomagus Mini Milling Drum Extra Fine	SBM-1002-34
Noviomagus Mini Milling Drum Extra Extra Fine	SBM-1002-35
Noviomagus Mini Milling Drum Extra Extra Extra Fine	SBM-1002-33
Sterilisation Tray Mini Mill	SBM-1002-36
Clamp Noviomagus Mini Mill	SBM-1002-37

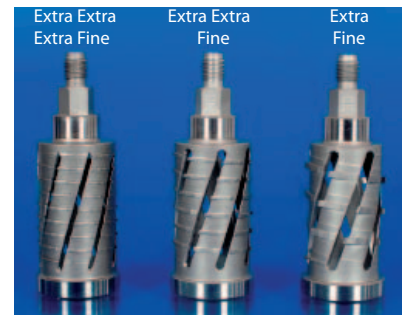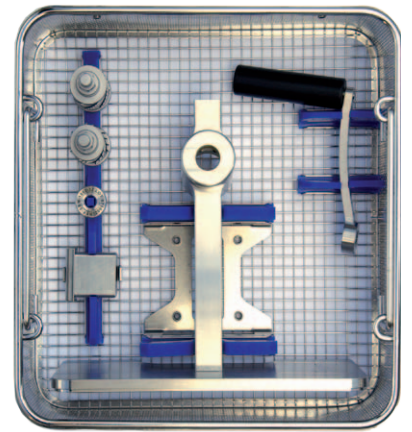

The Noviomagus Mini Mill is CE marked and FDA approved. US Patent No. 6,318,651. No rights can be derived from the information in this brochure. Copyright: Spierings Orthopaedics B.V.

Other Noviomagus products manufactured by Spierings Orthopaedics B.V. are: Noviomagus Bone Mill, Noviomagus Revision Meshes and Noviomagus Femoral Head Reamer Set.

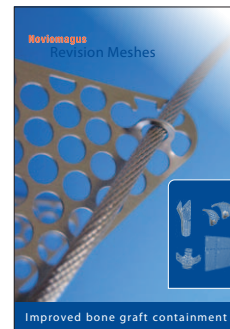

Manufactured by:
Spierings Orthopaedics B.V.
Madoerastraat 24
6524 LH Nijmegen
The Netherlands
T: +31 -(0)24-3501603
F: +31 -(0)24-3501604
E: info@spierings.biz
I: www.spierings.biz